



HIGH SPEED REFRIGERATED MICRO- CENTRIFUGE, FLOW R24

Please read the User Manual carefully before use, and follow all operating and safety instructions

User Manual
English

User Manual

EN

HIGH SPEED REFRIGERATED MICRO CENTRIFUGE, FLOW R24

Preface

Thank you for purchasing our product. Users should read this manual carefully, follow the instructions and procedures, and be aware of all preventive measures when using this instrument.

Service

If help is needed, you can always contact your dealer or Labbox via www.labbox.com.

Please provide the customer service representative with the following information:

- Serial number
- Description of the problem
- Your contact information

Warranty

This instrument is guaranteed to be free from defects in materials and workmanship under normal use and service for a period of 24 months from the date of invoice. The warranty is extended only to the original purchaser and shall not apply to any product or parts that have been damaged due to improper installation, improper connections, misuse, accidents, or abnormal conditions of operation.

For claims under the warranty, please contact your supplier.

CONTENTS

1. INTRODUCTION	4
2. INTENDED USE	4
3. SYMBOL	4
4. FEATURES	5
5. STANDARD ACCESSORIES	6
6. TECHNICAL SPECIFICATIONS	6
7. SAFETY PRECAUTIONS	7
8. INSTALLATION	10
9. STANDARD PARTS LISTING	11
10. USER INTERFACE AND DISPLAY	11
11. ROTOR INSTALLATION	13
12. OPERATING THE CENTRIFUGE	16
13. REMOTE OPERATION & PROGRAMMING	16
14. MAINTENANCE AND CLEANING	21
15. TROUBLESHOOTING	23
16. WARRANTY STATEMENT	26
17. PRODUCT DISPOSAL	27
18. TRANSPORTATION & STORAGE	27

1. INTRODUCTION

This manual provides important safety information for this centrifuge. It should be kept near the centrifuge for quick & easy reference. This centrifuge is equipped with maintenance free brushless DC motor and microprocessor controlled max speed up to 15,000 RPM and a safety feature of imbalance detection. This centrifuge is equipped with refrigeration which can be maintained as low as up to -20°C up to 40°C.






2. INTENDED USE


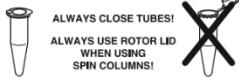

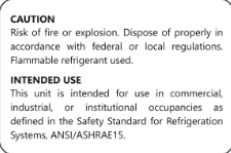
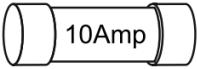


This centrifuge is designed to separate, sediment, spin down aqueous solutions & solvent suspensions of differing densities in compatible sample containers.

NOTE: Before using the centrifuge, please read this user manual carefully. This user manual is intended to assist with the operation and care of the unit & is not a document which aids in repair. For repair, please contact the supplier.

3. SYMBOL

Information on the device

Information	Meaning	Location
	Warning Observe the safety-relevant information in the operating manual	Keypad
	Read the operating manual.	Right side of the device
	Warning The device contains the flammable R290 Refrigerant.	Right side of the device
	Warning Risk of cryogenic burns.	Upper side of the device, under the centrifuge lid
	Warning Risk of hand injuries.	Upper side of the device, under the centrifuge lid

	Risk of damage to health when handling infectious liquids of pathogenic germs.	Upper side of the device, under the centrifuge lid
	Always seal the tubes. If you are inserting centrifuge tubes, use the rotor lid to seal the rotor	Upper side of the device, under the centrifuge lid
	Tighten the rotor nut using the rotor key provided.	Upper side of the device, under the centrifuge lid
	Fire hazard or explosion hazard. Dispose of the device in accordance with the applicable laws and regulation. The device contains a flammable refrigerant. Use the device in accordance with the intended use.	Right side of the device
	Fuse	Right side of the device (Near power socket)
	High Voltage	Right side of the device (Near power socket)
	Waste electrical	Right side of the device (In barcode sticker)

4. FEATURES

- Maximum volume of 44 x 1.5/2 ml
- Temperature ranges from -20°C to 40°C
- Temperature of 4°C maintained even at high speed
- Imbalance detector with auto cutoff for safe operation
- Brushless DC motor for maintenance free performance

- Large LCD display with speed and time indicator
- Speed setting from 500-15000 RPM
- One press changeover from RPM to RCF function
- Digital countdown Time
- Reduction adaptors for 0.1/0.2 & 0.5ml tubes
- Emergency lid release during power failure
- Automatic internal diagnosis & error display
- Silent operation

5. STANDARD ACCESSORIES

- T - Allen wrench
- Adaptors for 0.4/0.5 ml tubes
- Adaptors for 0.1/0.2 ml tubes
- Product user manual and warranty card
- Power cord
- Grease Tube
- Fuse 10A

NOTE: Rotors must be ordered separately. (optional accessory)

6. TECHNICAL SPECIFICATIONS

Motor Type	Brushless DC Motor
Rotor capacity	24 x 1.5/2ml tube
	44 x 1.5/2ml tube
	8 x 4 x 0.1/0.2ml PCR tube/strips
	8 x 5ml tube
Speed Setting	Variable 500 - 15000 RPM
	22388 x g (Max RCF)
Speed Step	100 RPM
Speed Accuracy	± 100 RPM
Temperature Accuracy	± 1°C

Temperature Range	-20°C to 40°C
Run Time	30s to 999min & 59s
Acceleration Time	< 72 seconds
Deceleration Time	< 74 seconds
Noise Level	< 60 dB
Ambient Temperature	≤80%
Permissible Relative Humidity	20 to 80%
Size (W x D x H)	325 x 720 x 305 mm
Weight	43 kgs
Refrigerant	R290
Total mass of Refrigerant	90 grams
Maximum allowable pressure (PSI) under normal condition	260 psi
High end pressure (PSI)	160 psi
Low end pressure (PSI)	10 psi
Altitude	Use up to altitude of 2000 m above MSL
Pollution Degree	2
Environment	For indoor use only
Power Supply	230 VAC, 50Hz
Power Consumption	650W
Safety Fuse Rating	10Amp

(* RCF varies for each rotor congregation - As mentioned in point no.11 (rotor installation).)

7. SAFETY PRECAUTIONS



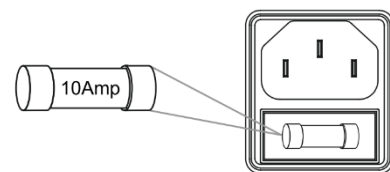
Read all safety & usage information provided in this manual carefully before using the device.

- Never use the centrifuge in any manner not specified in this manual.
- Using equipment in any manner not specified in this manual or by the manufacturer, will void the warranty.
- Always use recommended original rotors and spare parts for best result &

product safety.

- Never move the centrifuge while the rotor is spinning.
- The rotor and the rotor lid must always be securely fastened. If the centrifuge shakes unusual noise during operation, the rotor or rotor lid needs to be checked. Switch OFF the device immediately by pressing STOP, check rotor fit & fasten well.
- The rotors must be loaded symmetrically. Each tube should be counterbalanced by another tube of the same weight.
- Do not use the centrifuge or rotor that have not been correctly installed or damaged.
- Repairs must only be performed by authorized service technician. Using incorrect rotors & wrong spare parts will void the warranty.
- Centrifuge may be used for the specified applications only. It must not be operated in hazardous or flammable environment and must not be used to centrifuge explosive or highly reactive substances. Also do not place the potential hazardous materials within the clearance area/envelope.
- If liquids are spilled on the rotor or rotor chamber, the centrifuge must be cleaned carefully and properly before being used again.
- Prior to centrifugation, the tubes should be visually inspected for material damage. Damaged tubes may not be centrifuged. This is because broken tubes can, in addition to sample loss, create imbalance which can result in further damage to the centrifuge and accessories.
- The maximum capacity of tubes must not be exceeded corresponding rotor capacity. Do not use liquid with density higher than 1.2g/ml for full load operation.
- Do not lean on the equipment. It may damage the equipment or harm the operator.
- When moving the centrifuge from a cold room to a normal room, run the centrifuge for 30 minutes beforehand in the cold room to avoid condensation. Alternatively, allow it to warm up in the lab for at least 3 hours before use, but do not plug in the centrifuge in order to prevent damage by condensation.

- Be sure to close the tubes lid tightly prior to centrifugation. Open tubes lid can be torn off during centrifugation and can damage the rotor lid or centrifuge.
- Rotors and rotor lids are high-graded components which are subject to extreme mechanical strain. Even slight scratches or tears can lead to serious material damage. Ensure to check the rotor for any signs of damage before use. Rotor & rotor lid showing visible signs of corrosion or mechanical damage should not be used, contact your local supplier or manufacturer.
- Do not fill tubes while they are in the rotor. Liquid spillage may harm the device.
- In the event of contamination caused by aggressive agents, the rotor must be cleaned immediately using 70% IPA. This is particularly important for the bores of the tubes. If any damage is seen, contact the service technician.
- Before using cleaning or decontamination methods other than those mentioned by the manufacturer, contact the manufacturer to ensure that the intended method will not damage the centrifuge.
- The power cord given with centrifuge unit is designed to use for that particular centrifuge. Do not use any other power cord, using any other power cord/adaptor may damage the centrifuge and will void the warranty.
- Safety Fuse is provided of 10Amp configuration, which can be replace by the operator. The same will protect the machine circuit during an overload fault or overload.
- For safety we have provided protective earthing with power supply. Make sure power supply is earthed.
- Be sure to close the tubes lid tightly prior to centrifugation.
- For safety and air-tight sealing to minimize the air noise, gentle force needs to be applied on the lid to lock properly.



WARNING: Keep all ventilation openings in the enclosure, or in the structure for building-in, clear of obstruction.

WARNING: Do not use mechanical devices or other means to accelerate the defrosting process, other than those recommended by the manufacturer.

WARNING: Do not damage the refrigerant circuit.

- The R290 refrigerant used in the device is flammable. If the refrigeration cycle is defective, the refrigerant may escape and form explosive mixtures with the ambient air.
- Observe the regulations that apply to your laboratory.
- Make sure there is a sufficient volume of air at the location.
- Ensure adequate ventilation of the device at the location.

8. INSTALLATION

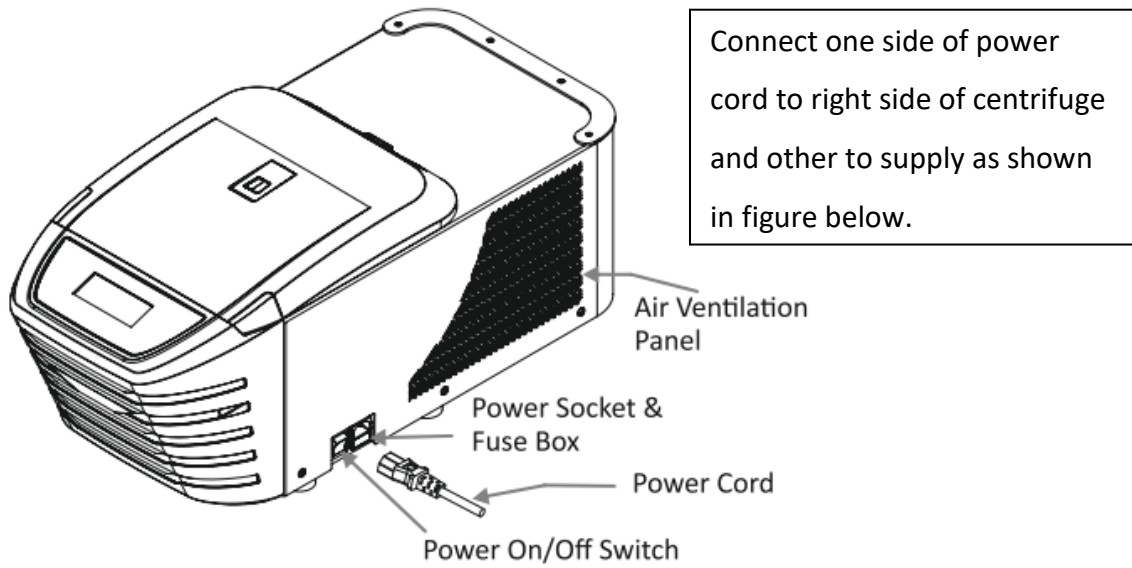
Open the box, then remove the packaging and gently take the centrifuge out of the box. Before using this centrifuge open centrifuge & remove all packaging from the rotor chamber & ensure rotor is firmly tightened. The user manual and accessories should be kept near the centrifuge. Please keep all packaging in safe storage for warranty purposes.

MOUNTING


Place the centrifuge on a flat and leveled surface, ensure that the four feet of this centrifuge stand on the surface firmly. Avoid installing on a slippery surface or surface prone to vibration.

1. Ideal ambient temperature is $20^{\circ}\text{C} \pm 5^{\circ}\text{C}$, avoid placing the centrifuge in direct sunlight.
2. Keep clearance of at least 50 cm on both sides and at least 30 cm behind it to guarantee cooling efficiency.
3. Keep away from heat or water to avoid sample temperature issues or centrifuge failures.
4. Do not place the centrifuge in any area where operating the unit may be difficult.
5. Ensure there are no source of ignition (flame or sparking electrical components) within 3 meters of working area.


CONNECTING POWER CORD



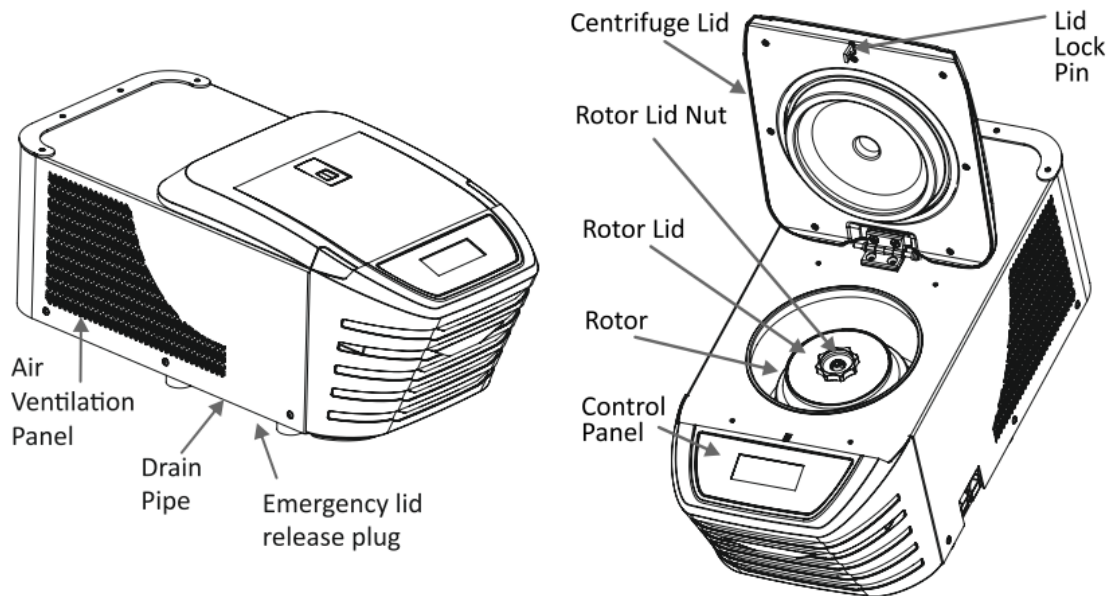
9. STANDARD PARTS LISTING



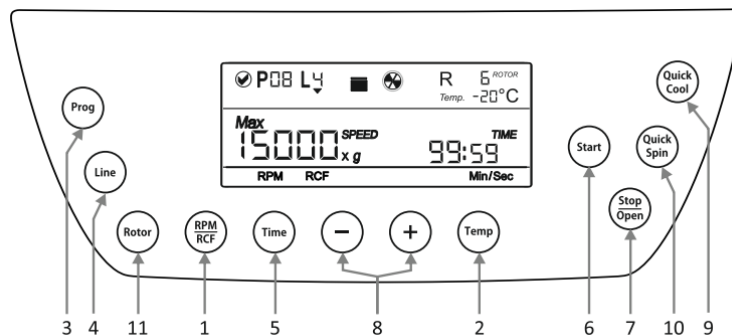10. USER INTERFACE & DISPLAY



Item	Name	Function
1	RPM/RCF	Press button to select SPEED value in rpm mode. And press +/- key to set desired run RPM. Press and Hold to switch to RCF mode. Press +/- key to set desired RCF.
2	TEMP	Press button to set the temperature value in a range from -20°C to 40°C. And press +/- key to set desired temperature.
3	PROG	Press button to enter in Program Mode and to select any Program out of the 99 available programs.
4	LINE	Press button to select line of any specific program out of 4 available lines in each program.
5	TIME	Press button to select TIME mode. And press +/- key to set desired run TIME.
6	START	Press button START to start the centrifuge.
7	STOP/OPEN	Single press Stop/Open button will stop motor. Second press Stop/Open button will open the lid once the motor stop. Also used for lid open if centrifuge motor is not running.
8	+/- Keys	Use to set speed, Temperature and time values. Press "+" to increase and "-" to decrease values.
9	QUICK COOL	To start a temperature control run at optimized (5000) RPM to achieve maximum cooling.
10	QUICK SPIN	To set short spin on a set speed. Press short spin continues press for Quick spin operation.
11	ROTOR	To select among multiple rotor listed below.

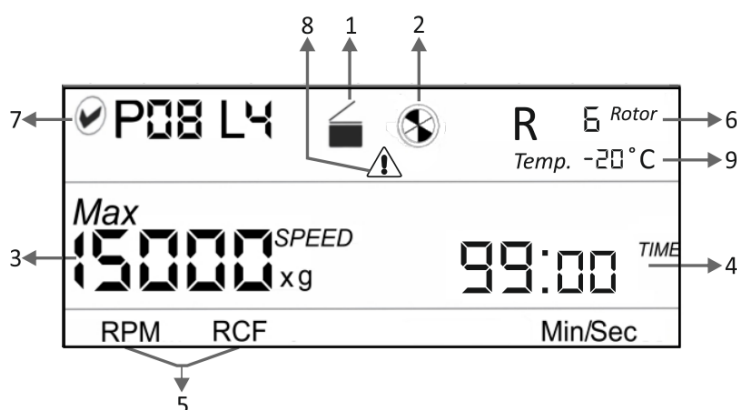

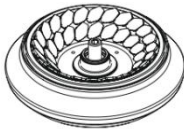



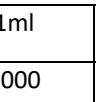

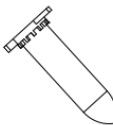
Item	Symbol	Description
1		Indicates lid status. Left image = lid close and Right image = lid open
2		Indicates centrifuge status. When centrifuge is running the symbol rotates and when centrifuge is not running the symbol is static.
3		RPM indicates the speed value at which centrifuge is running. Xg indicates the value in RCF mode.
4		Countdown timer. Indicates the time for which the centrifuge will run. Indicates the time in Min/Sec mode.
5		Indicates rpm or RCF mode and shows corresponding values.
6		Indicates selected rotor.
7		Indicates selected program number and the line number specific to the program.
8		Warning symbol. It appears when some error occurs.
9		Chamber temperature indication.

11. ROTOR INSTALLATION

ROTORS AND ACCESSORIES

The table below shows rotors compatible with the centrifuge and the Max. RCF with different tubes and adaptors.

24 x 2ml rotor (standard supply)	1,5/2.0ml tubes	0,5 ml Adaptor	0,1/0,2ml Adaptor		PCR strip rotor (Optional)	0,2ml PCR strip
Max. Speed	15000	15000	15000	15000	Max. Speed	15000
Max. RCF x g	21420	18842	16103	16930	Max. RCF x g	18461

44 x 2ml rotor (optional)	1,5/2.0ml tubes	0,5 ml Adaptor	0,1/0,2ml Adaptor		8 x 5ml rotor (Optional)	5ml tube
						
			0,1ml	0,2ml		
Max. Speed	15000	15000	15000	15000	Max. Speed	15000
Max. RCF x g	22388	19810	17070	17898	Max. RCF x g	18461

Note: Always use recommended accessories for best results and product safety.

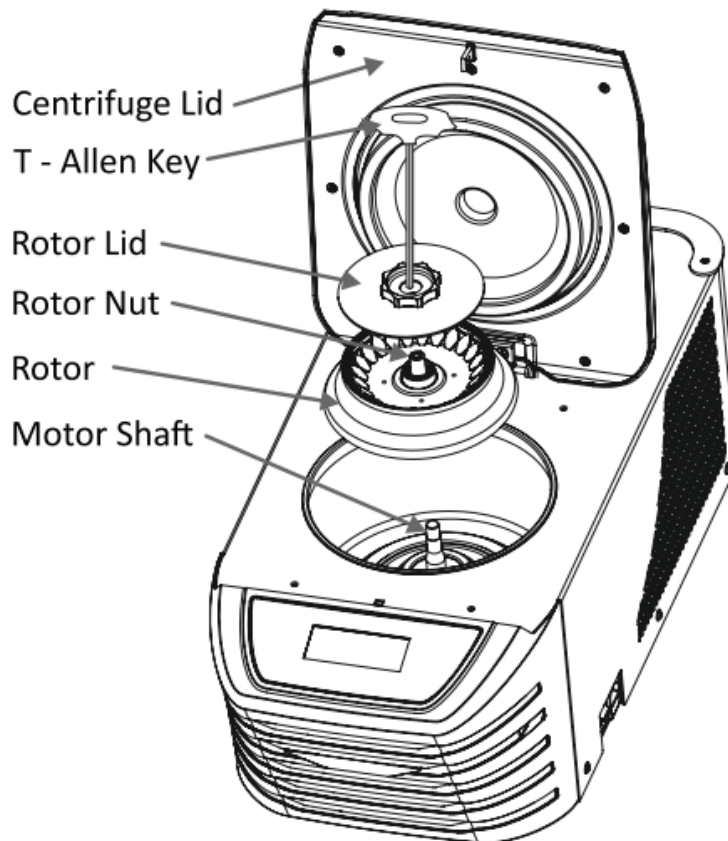
Optional rotors & additional accessories must be ordered separately.

Note: Rotor maximum service file 50.000 cycle (from first use).

REMOVING THE ROTOR

1. Do not remove or loosen the rotor lid before attempting to remove the rotor.
2. Using the T-Allen Key, loosen the rotor nut by turning counterclockwise. Do not try to pull the rotor, the rotor will come up automatically.
3. Once the rotor nut is loosened completely, pull up the rotor vertically.

REPLACING THE ROTOR




1. To replace or install the rotor, take the rotor and load vertically onto the shaft.
2. Place the rotor nut in the center hole of the rotor onto the motor shaft.
3. Put T-Allen Key into the rotor nut and turn clockwise to tighten and counterclockwise to loosen the rotor.
4. After properly fastening the rotor, place the rotor lid on the rotor lid nut by hand and rotate the rotor lid clockwise.


NOTE:

- 1) Check the rotor is firmly tightened before running the next program.
- 2) Do not remove or loosen the rotor lid before removing the rotor.

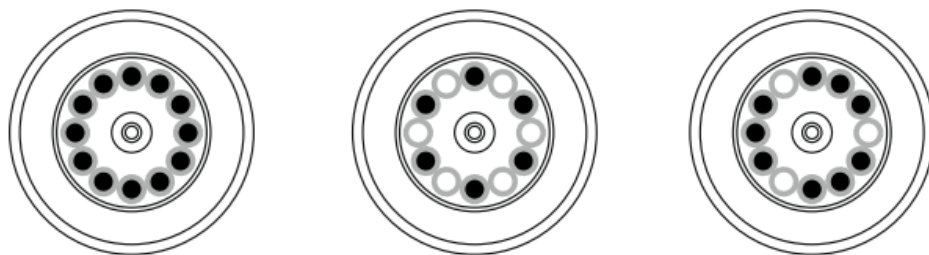BALANCING THE ROTOR




1. Always balance the rotor before centrifugation. Above are examples of properly balanced rotors.
2. The samples in the tubes should be of equal volume.
3. If the tubes are not loaded correctly - vibration or imbalance can occur, which can cause serious damage to the centrifuge.
4. If the tubes are not loaded symmetrically then the imbalance detector will cut off the running centrifuge for device & user safety. This will stop the centrifuge and Err 55 will be seen indicating tubes are not loaded symmetrically. To resume operation, load tubes symmetrically & restart the centrifuge.
5. Incorrect method of loading tubes in centrifuge rotor:

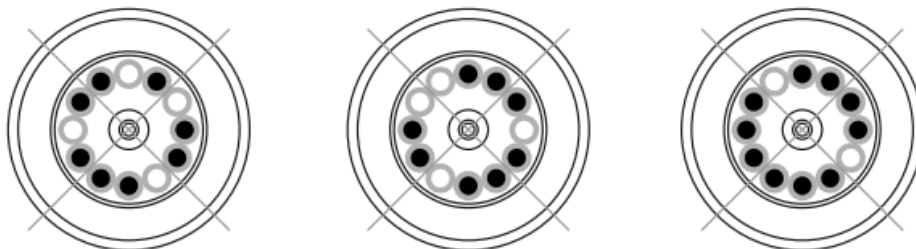

12. OPERATING THE CENTRIFUGE

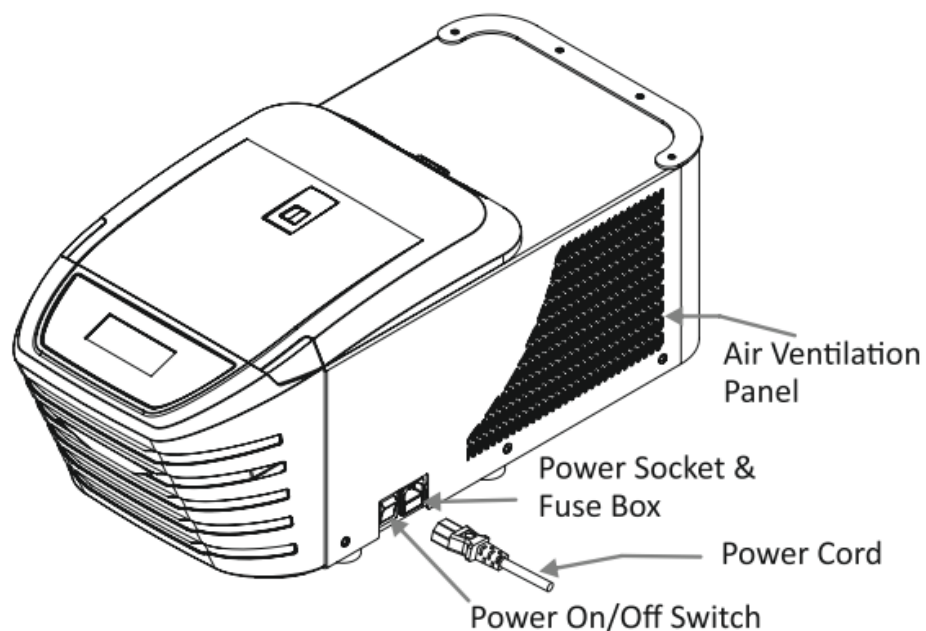STARTING THE CENTRIFUGE



After connecting the power cord, switch ON the power at the right side of the centrifuge. Make sure to check the rotor fit before use. Centrifuge will not operate with an open lid.
NOTE: Maintain a gap of 3 seconds between turning off and turning on the centrifuge. DO NOT turn off and turn on the centrifuge instantly.

13. REMOTE OPERATION & PROGRAMMING

SWITCH ON THE CENTRIFUGE

After connecting the power cord. Switch ON the main power supply & then Switch on the power switch located on the right side of the instrument. Make sure to check the rotor fitment before use. Centrifuge will not operate with open lid.

NOTE: Maintain a gap of 3 seconds between switch OFF and switch ON again. DO NOT switch OFF and ON again instantly.

SPEED SETTING

After closing the centrifuge lid, press “RPM / RCF” to select speed setting in RPM / RCF mode. Press “+” to increase the speed value and press “-” to decrease the speed value. The minimum and maximum speed (RPM) of the centrifuge is 500 rpm and 15000 rpm respectively. You can quickly change the increment / decrement value by rotating setting

knob quickly. Value is accepted when left idle indicated by multiple blinks of the set value.

The speed can also be changed while the centrifuge is under operation. Press the RPM / RCF button & use setting knob to change speed as instructed above. Changing the speed between the ongoing centrifugation will run the centrifuge at updated speed for the remaining time indicated by the timer.

TIMER SETTING

Press “TIME” button to set the run time. After selecting press “+” to increase the timer value and press “-” to decrease the timer value. Value is accepted when left idle indicated by multiple blinks of the set value. The centrifuge timer can be set from 30s to 999mins and 59s. The timer in the centrifuge is countdown timer and time in the display will be in “Min/Sec”. The same will be shown on the display.

1. To set the timer in minutes, pressing the time button “ONCE”. For example: if time is 05:00 (5 Mins 0 Sec), then next time will be 06:00 (6 Mins 0 Sec).
2. To set the timer in seconds, pressing the time button “TWICE”. For example: if time is 05:00 (5 Mins 0 Sec), then next time will be 05:01 (5 Mins 1 Sec).
3. Pressing time button more than 2 times will start the process from point 1 again.
4. The input will be accepted if we leave the setting knob idle for 3 seconds. The value be saved automatically after five blinks.

TEMPERATURE SETTING

Press temperature button for setting the temperature from -20°C up to 40°C. The parameters can be set when the Temperature value blinks by pressing “+” key for increasing the temperature & “-” to reduce the temperature. The values will automatically save after screen stops blinking.

QUICK COOL SETTING

PRESS Quick Cool Key, indicated by (cubic rotation) on top right (Indication will be in dotted rotation in place of rotor indication). During quick cool operation rpm: - FIX 5000 & TEMP: -On set value. Press START button for centrifuge to start at FIX RPM & set TEMP. After reaching set TEMP, motor stops thus RPM becomes zero, but compressor continues to keep the chamber temperature at the set value. The quick cool works for 1 hour only. (for indication of refrigeration running process, TEMP display will blink) If

timers stops and TEMP does not reach to set value, TIME will switch-over from beginning. "RPM/RCF, Timer, Rotor, Line, Prog", "Quick Spin" button does not work during "Quick Cool" operation.

SHORT SPIN OPERATION

Short Spin Centrifugation is the feature for short/pulse/quick run. It will run as long as the button is pressed. Set rotational speed prior to short spin as required. During short spin the timer will be in incremental mode. After releasing short spin button, the time in the display will show duration of short spin.

ROTOR SELECTION


Different rotors have different radius which means that RCF for that rotor would also be different. Proper selection of rotor is mandatory to achieve proper RCF. Press "Rotor" button then rotate setting knob to select the rotor from available 4 rotors. Press "Rotor" button. See below chart to identify rotor no. from its volume:

Rotors			
Rotor No.	Description	rpm	RCF at max rpm
Included	24 slot x 2ml rotor	15000	21420 g
CCAN-A11-001	8 slot x 5ml rotor	15000	18461 g
CCAN-A12-001	44 slot x 2ml rotor	15000	22388 g
CCAN-A13-001	4 x 8 PCR rotor	15000	15397 g

Press "+" & "-" button to change the value/no of rotor to set that rotor no or value will be stored automatically if no button is pressed after 6 seconds of the adjustment. Then set max value of speed and acc/dcc to be changed as per the rotor selection. (Max. Speed (RPM/RCF) and Acc/Dcc value for all the rotors have been added in the software) Now set the desired speed (RPM/RCF), time & temp and press START button to store and run the selection.

START AND STOP OPERATION

Press "START BUTTON" to start operation and press "STOP/OPEN BUTTON" to stop the ongoing operation. When the centrifuge is running the symbol "⚙️" will rotate. Pressing the "STOP/OPEN BUTTON" will stop the operation. Once operation is stopped, press Button again to open the centrifuge lid. If the time gets over, centrifuge will stop

automatically. When the centrifuge is not running the symbol “” will be idle. To open the lid in non-operating stage, press the “STOP/OPEN BUTTON”.

After completing run, before any another operation, it is mandatory to open the lid once before starting new operation. Centrifuge will not do 2nd operation if lid is not opened and closed for at least one time after completion of 1st operation.

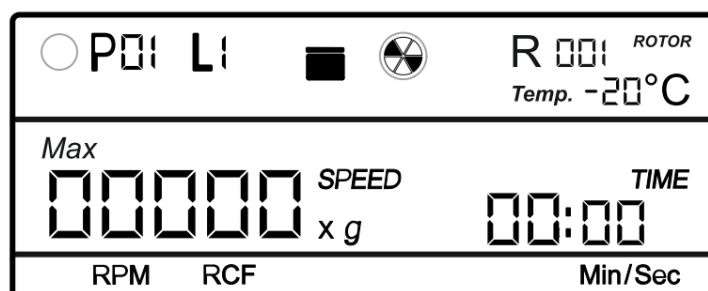
NOTE: It is mandatory to open the lid once after completion of operation for 2nd operation. Centrifuge will not start if lid is not opened and closed after completion of operation.

PROGRAM MODE


Long press “PROG” button to enter program mode. Short press select program and rotate setting knob (clockwise & anticlockwise) to select desire program out of 99 available programs. User can select any program from 1 to 99 programs and can have user specific setting parameters for all programs. Each program has 4 lines indicated L1, L2, L3, L4. Press “LINE” button to enter in line selection mode and rotate setting knob to select one out of 4 lines.

Every line of the program can have different values for all setting parameters. User can set and save different speed, different time, different temperature for any line of program. Only Rotor will not remain save for any program. Rotor selection is not Program specific. Rotor selection is universal and it should be set before setting any other parameters.

On first time usage, all programs will have zero (0) values and once the “PROG” button is pressed below display will appear indicating selection of program mode.




As per above image, the circle o left of “P01 L1” is blank indicating no values in any line of specific program. Once any line of Program is filled or set, below change of display will appear.

The circle  left of “P01 L3” will have a tick mark indicating that any one line of the


specific program is filled. If any line leading to selected line is set, then the arrow below the line number will appear.

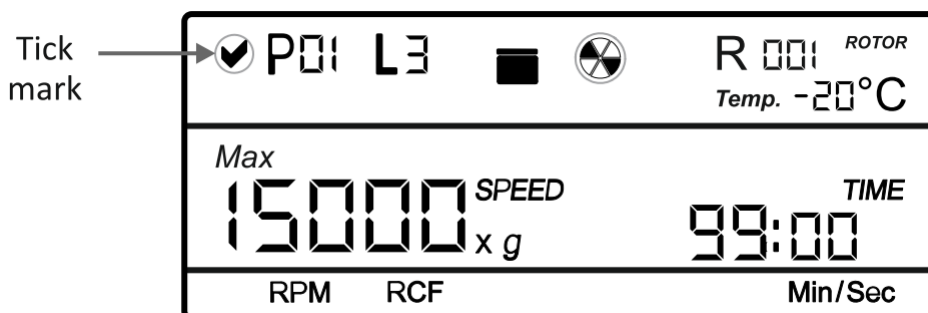If no leading line is filled or last line L4 is selected, then the arrow below Line will not appear. It indicates there is no further operation left as leading lines of the specific program are empty.



Program and selected rotor remain active for all programs until new rotor number is not selected. Select the rotor, Press “PROG” button to enter in Program mode, select specific program using setting knob, Press “LINE” button to select specific line and set different parameter for that particular line of program. Values for any parameter gets save after 3 blinks. Once Program is set, press “START” button to start program.

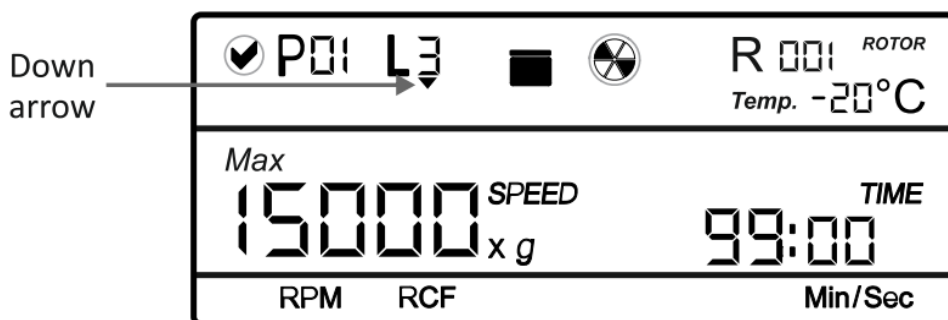After completion of any operation, it is mandatory to open the lid at least once for another operation. Open the lid, close it again and press START for another operation.



NOTE:

- 1) Rotor is not program specific. Selected rotor will be active for any new operation.
- 2) It is mandatory to open the lid once after completion of operation for 2nd operation. Centrifuge will not start if lid is not opened and closed after completion of ongoing operation.

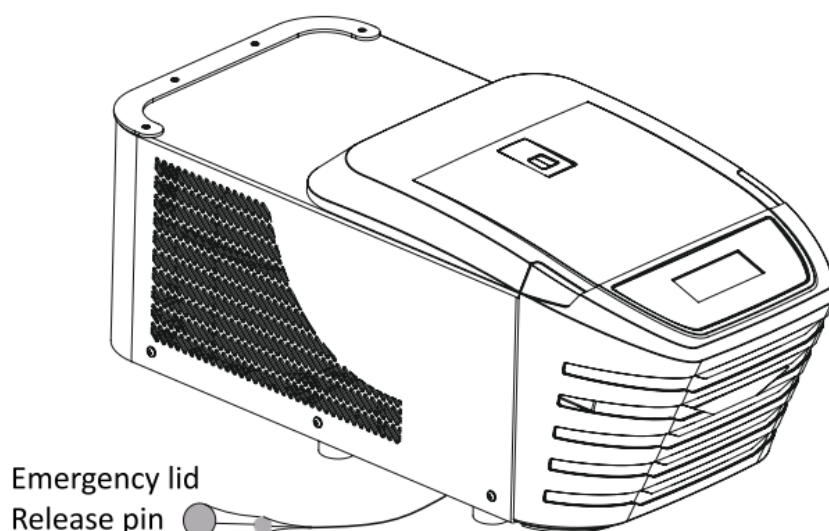
IMBALANCE DETECTION

The centrifuge is equipped with an imbalance detection safety feature. When the rotor is not loaded symmetrically, the imbalance detector gets activated and will cut off the centrifugation. The error “Err 55” will be shown on the display. First correct the imbalance load then switch OFF & switch it ON again. The values will be same as set before imbalance. The imbalance detection feature cannot be deactivated, as it is factory fitted.

OPENING CENTRIFUGE LID IN POWER FAILURE

Disconnect the centrifuge from the main power supply. Wait until the rotor has come to a standstill (this may take longer time). Once the rotor has stopped, Pull the Emergency Lid release Thread located at the bottom left (to your right) of the machine. This will open the centrifuge Lid.

NOTE: This method of opening the lid should be used only in case of Emergency or power failure.



14. MAINTENANCE AND CLEANING

- The rotor and the outside of the centrifuge should be cleaned regularly with a mild wet (with water) cloth.
- Ensure that while cleaning the unit is not plugging in.
- Wear protective glove & safety glass while operating & cleaning the device.
- The brushless motor in the centrifuge requires no routine maintenance. Any required service should be performed by authorized, qualified personnel only.

Repairs performed by unauthorized personnel may void the warranty.

- Always keep the centrifuge housing, rotor chamber and rotor clean. All parts should be wiped down periodically with a soft cloth. For more thorough cleaning, use a neutral cleaning agent (Ph between 6 and 8) and clean with a soft cloth. Excessive amounts of liquid should be avoided.
- **NOTE:** liquid should not come in contact with the motor.
- After cleaning, ensure that all parts are dry before re-use.
- Regularly cleaning of the rotor is important.
- If the rotor chamber needs cleaning, clean with cloth or sponge moistened with a neutral detergent solution.
- Do not place the rotor into the cleaning solution.
- If corrosive, toxic or pathogenic bacteria are accidentally spilled in the rotor or rotor chamber the centrifuge must be decontaminated thoroughly.
- Grease tube has been provided as part of standard accessory. Minor Grease has to be applied on the Motor shaft (Threaded Parts, Grooves) at every 2 weeks interval to prevent rotor & shaft jam. This also aids ease installation process in case of regular rotor change.
- Switch off the device and disconnect it from the mains/power line before starting cleaning or disinfection.
- Do not connect the device to the mains/power line unless both the inside and outside of the device are completely dry.

Checking the device:

- Check the device, main/power cord and motor shaft for visible damage. If the device is damaged, take it out of operation. Inform the authorized service technician.
- Check the rotor chamber for corrosion, In case of corrosion in the rotor chamber, inform the authorized technician.

Checking the rotor and accessories:

- Check the rotors and accessories for damage. Take damaged rotors and damaged accessories out of operation.
- Check the service life of rotors and accessories. Replace rotors and accessories that have exceeded their service life.

Documenting the service life:

- The service of rotors, rotor lids and accessories may be limited. Information on how long you may use rotors can be found in the mentioned technical data.

Document the following data for rotors, rotor lids and accessories:

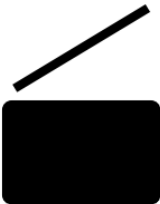
- Data of initial setup.
- Number of centrifugation cycles.
- Number of autoclaving cycles.

Fuse Replacement Instruction:

- Disconnect the device from the electrical system.
- Remove the fuse from power socket. In some cases, you may need a small minus screwdriver to unscrew the fuse holder cap.
- Look at the fuse wire. If there is a visible gap in the wire or a dark or metallic smear inside the glass, then the fuse is blown and needs to be replaced.
- Replace the fuse with one that is the same. Make sure to note the fuse ampere and voltage ratings.

15. TROUBLESHOOTING

This centrifuge has a self – diagnostic function. If a problem occurs, an error/warning code will be displayed on the display screen, and the operator can determine the malfunction with the warning code below.

ERROR	PROBLEM	SOLUTION
No display	No main power connection.	Power check & proper plug-in of main cable at both ends.
	Power failure.	Check the mains fuse.
	Improper connection.	Connect adaptor properly.
	Fuse Blown of unit.	Replace the Fuse of unit.
	Lid not closed correctly.	Close lid correctly and restart unit.
	Error with lid closing and opening mechanism.	Contact service.
Error 1	Latch jammed or any limit switch of latch got damaged.	Contact service representative.
Error 41	Selected temp value not reachable for selected set speed.	Set temp value will be taken automatically that it can be archived after 1 hour. Error indication will be displayed every 15 second. restart unit.
Error 42	Temperature is out of control.	Error indication will be displayed every 15 second. If persist in every operation consistently, contact service representative. restart unit.
Error 43	Over temperature inside chamber.	Turn off the centrifuge and wait until temperature goes down. restart unit.
Error 44	Temperature sensor failure.	Contact service representative.

Error 52	Rotor stuck or incorrect operating voltage.	Turn off the centrifuge, check rotor fitment or apply correct 230 VAC \pm 10VAC operating voltage. restart unit.
Error 55	Rotor not loaded symmetrically.	Load rotor symmetrically & restart centrifuge. restart unit.
Centrifuge lid cannot be opened.	Rotor is still spinning.	Wait for the rotor to come to a stop.
	Power Failure.	Emergency lid release after rotor stops.
Centrifuge shakes during acceleration & exceptional running noise.	Rotor not loaded symmetrically.	Load rotor symmetrically & restart operation.
	Either a broken tube, damage to the rotor or motor is cause for run noise.	Replace broken tube. For damaged rotor/motor contact service representative.
	Rotor damaged.	Remove & change rotor.
Display Error.	Loose connection of display.	Contact service representative.
Power tripping.	Cable not fit properly.	Remove cable and connect properly.
System gets hang.	Electronic error.	Switch OFF centrifuge and then switch it ON again. If the error still shows, contact service representative.

IMPORTANT NOTE

1. Maintain a 3 second gap between restarting the centrifuge. Instant ON-OFF can lead to a reset, erasing last run memory.
2. If the motor gets hot due to which there will be a fluctuation in speed value, then allow the centrifuge to cool for at least 30 minutes. Do not do any operation for 30 minutes.
3. Do not use liquids with density higher than 1.2g/ml for full load operation.

16. WARRANTY STATEMENT

This product is warranted to be free from defects in material and workmanship. Your product will be duly repaired upon prompt notification in compliance with the following conditions:

This warranty is valid only if the product is used for its intended purpose and within the guidelines specified in this instruction manual. This warranty does not cover damage caused by accident, neglect, misuse, improper service, natural forces or other causes not arising from defects in original material or workmanship. This warranty does not cover any incidental or consequential damages, commercial loss or any other damages from the use of this product.

The warranty is invalidated by any non-factory modification, which will immediately terminate all liabilities on us for the products or damages caused by its use. The customer shall be responsible for the product or use of products as well as any supervision required for safety. If requested the products must be returned to the distributor in well packed and insured manner and all shipping charges must be paid.

NOTE: Some states do not allow limitation on the length of implied warranties or the exclusion or limitation of incidental or consequential damages. This warranty gives you specific legal rights. This warranty is given expressly in lieu of all other warranties, expressed or implied.

The purchaser agrees that there is no warranty of merchantability or of fitness for any intended purpose and there are no other remedies or warranties, implied, which extend beyond the description on the face of the agreement. This warranty is only applicable to the original purchaser.

Products received without proper authorization will not be processed for warranty or service. All items returned for service should be sent with postage prepaid in the original packaging or other suitable packaging, padded to avoid damage. We will not be responsible for damage incurred by improper packaging.

For details regarding the warranty period, please refer to the warranty card provided with the product.

NOTE: This warranty is valid only if the warranty is registered with the supplier within 30 days from the date of purchase.

For your reference, make a note of serial number, date of purchase and supplier here.	
Serial No.:	Purchase Date:
Supplier:	

17. PRODUCT DISPOSAL

In case the product is to be disposed of, the relevant legal regulations are to be observed. Information on the disposal of electrical and electronic devices in the European Community.

The disposal of electrical devices is regulated within the European Community by national regulations based on EU Directive 2012/19/EU on waste electrical and electronic equipment (WEEE). According to these regulations, any devices supplied after 13.06.05 in the business-to-business sphere, to which this product is assigned, may no longer be disposed of in municipal or domestic waste. They are marked with the following symbol to indicate this.

As disposal regulations within the EU may vary from country to country, please contact your supplier if necessary.

Caution! Fire hazard or explosion hazard

The device contains a flammable refrigerant. If the refrigerants used come into contact with sparks, they will ignite. Persons may be injured.



- Dispose of the device properly in accordance with national or local regulations.



18. TRANSPORTATION & STORAGE

This device is heavy & weighs approximately 43kgs. Care must be taken while lifting. Always take few people help to lift this instrument. The refrigerated centrifuge must be kept always in upright position while storing/transporting from place to place.

- Use only original packaging during transportation.
- For longer distance take transportation aid like hard trucks.
- Avoid knocking, harsh shaking or jolting the device.
- Always retain the packaging material & transportation protections for longer storage or transportation.

- The transportation conditions for the instrument are -25°C to 60°C with a Relative Humidity up to 80% & max pressure of 106kPa for both general & air transportation.
- For storage the instrument is recommended to be stored in original package. the instrument is -5°C to 45°C with a Relative Humidity up to 70% & max pressure of 106kPa.

Nota importante para los aparatos electrónicos vendidos en España

Instrucciones sobre la protección del medio ambiente y la eliminación de aparatos electrónicos:



Los aparatos eléctricos y electrónicos marcados con este símbolo no pueden ser eliminados en forma de residuos urbanos.

De conformidad con la Directiva 2012/19/UE, los usuarios de la Unión Europea de aparatos eléctricos y electrónicos, tienen la posibilidad de devolver sus RAEE para su eliminación al distribuidor o fabricante del equipo después de la compra de uno nuevo. La eliminación ilegal de aparatos eléctricos y electrónicos es castigada con multa administrativa.

Remarque importante pour les appareils électroniques vendus en France

Informations sur la protection du milieu environnemental et élimination des déchets électroniques :



Les appareils électriques et électroniques portant ce symbole ne peuvent pas être jetés dans les décharges.

En réponse à la réglementation, Labbox remplit ses obligations relatives à la fin de vie des équipements électriques de laboratoire qu'il met sur le marché en finançant la filière de recyclage de ecosystem dédiée aux DEEE Pro qui les reprend gratuitement (plus d'informations sur www.ecosystem.eco). L'élimination illégale d'appareils électriques et électroniques est punie d'amende administrative.

Nota importante per le apparecchiature elettroniche vendute in Italia

Istruzioni sulla protezione ambientale e sullo smaltimento dei dispositivi elettronici:



Le apparecchiature elettriche ed elettroniche contrassegnate con questo simbolo non possono essere smaltite come rifiuti urbani.

In conformità con la Direttiva 2012/19 / UE, gli utenti dell'Unione Europea di apparecchiature elettriche ed elettroniche hanno la possibilità di restituire i propri RAEE per lo smaltimento al distributore o al produttore di apparecchiature dopo averne acquistato uno nuovo. La rimozione illegale di apparecchiature elettriche ed elettroniche è punibile con una sanzione amministrativa.

www.labbox.com